

Bio-bin® Waste Disposal Containers



Bio-bin® Waste Disposal Containers

The ideal solution for biohazard or non-hazardous waste
REDUCE Cost, Weight, Transport, CO2 Emissions, Storage Space

- Space efficient, light weight fold-flat design
- Built-in lid and handles
- Temporary & permanent seals
- Stackable post use. UN 3291 Approved
- Printed user instructions & user ID/date information on lid
- Waterproof construction
- Built-in super absorbent Bio-matt™ to soak up any excess liquid in Bio-bin
- Autoclavable. Use with #202065 Autoclave Sterility Tape.
- Supplied with 2 x 2" Biohazard Waste Labels
- Reorder information printed on the Bio-bin

Visit
thomassci.com
 to view the
 Bio-bin
 video!

Catalog Number	Description	UOM	Price
1217C52	Bio-Bin, Loop Model, 2.5L	50/CS	\$146.69
1217C53	Bio-Bin, Pipette Model, 6L	40/CS	\$192.70
1217C54	Bio-Bin, Floor Model, 30L	10/CS	\$143.60
1217C55	Bio-Bin, U-Stand, for 6L	1 EA	\$30.83
1217C56	Bio-Bin, Z-Stand, for 6L	1 EA	\$30.83
1205X22	Bio-Bin, M-Stand, for 30L	1 EA	\$105.98

Bio-bins are a new cost effective and convenient alternative to rigid plastic waste collection and disposal containers. Ideal for use in academic, clinical, government and industrial laboratories for a wide range of non-sharps biomedical and laboratory waste that require special heat treatment such as autoclaving or incineration before disposal. Items such as microbiological loops, pipettes, tips, test tubes, swabs and vials are ideal for disposal by this method.

Thomas
 Scientific

800.345.2100 • 800.345.5232 fax
 PO Box 99 • Swedesboro, NJ 08085
www.thomassci.com • value@thomassci.com



Fold-flat

2.5L Bio-bin® is delivered folded flat

**Temporary Seal**

When lid is folded down press here for temporary seal

Permanent Seal

Peel off release paper and fold down lid to seal permanently

Printed user instructions

Simple instructions on how to seal the Bio-bin®

Built in lid

Push inwards here to close lid

Light weight

Weights up to 65% less than plastic alternatives

Biohazard

Supplied with 2" x 2" Biohazard labels for user application if necessary

**No sharps**

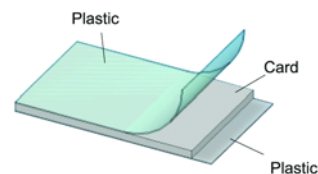
Only used for non-sharp infectious waste.

Stackable pre and post use

6L Bio-bin® is delivered stacked inside each other making them very space efficient. Once sealed down they can be stacked on top of each other.

**Water proof construction**

96% paper construction with a very thin coating of plastic laminated to make it water proof inside and out. Heat sealed joints prevent any leakage.

**Absorbent Bio-matt™**

Every Bio-bin® will come with a removable Bio-matt™. This will absorb liquid and turn it into a gel in minutes.



Thomas
Scientific