

SERF2 Polyclonal Antibody

Catalog number: 11691-1-AP


Size: 23 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:


SERF2; 4F5REL, FAM2C,
FLJ20431, FLJ37527, FLJ38557,
H4F5rel, HsT17089, Protein 4F5
related, SERF2, small EDRK rich
factor 2

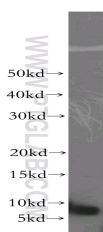

HepG2 cells were subjected to SDS
PAGE followed by western blot with
11691-1-AP (SERF2 antibody) at
dilution of 1:500

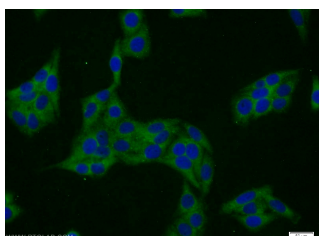

Immunofluorescent analysis of
HepG2 cells using 11691-1-AP (SERF2
Antibody) at dilution of 1:25 and
Alexa Fluor 488-conjugated
AffiniPure Goat Anti-Rabbit IgG(H+L)

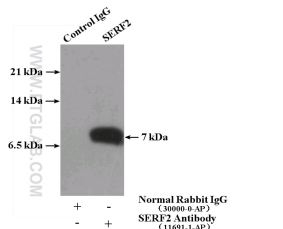

IP Result of anti-SERF2 (IP:11691-1-
AP, 4ug; Detection:11691-1-AP 1:500)
with HepG2 cells lysate 4000ug.

Background

SERF2 Belongs to the SERF family. The antibody is specific to SERF2.

Applications

Tested applications:	ELISA, WB, IF, IP
Species specificity:	Human; other species not tested.
Calculated SERF2 MW:	59aa, 7 kDa
Observed SERF2 MW:	7 kDa
Positive WB detected in	HepG2 cells, human adrenal gland tissue
Positive IP detected in	HepG2 cells
Positive IF detected in	HepG2 cells
Recommended dilution:	WB: 1:500-1:5000 IP: 1:200-1:2000 IF: 1:10-1:100

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag2305
GenBank accession number:	BC022326
Gene ID (NCBI):	10169
Full name:	Small EDRK-rich factor 2

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.