

DASH FLEX 24

PROGRAMMABLE STAT CENTRIFUGE

Part # 00-184-009-000



Always know when your STAT samples are ready with Drucker's DASH Flex 24 centrifuge. The unique flashing lid lights (patent pending) make it easy to track cycle status and keep TAT down. Customize settings and check every detail on the digital display, and save up to 10 custom cycles to produce high-quality STAT samples.

Easy to Monitor

LED lid lighting indicates when the centrifuge is ready (off), running (on), or done (flashing), reducing forgotten samples and keeping your average TAT down. The DASH Flex 24's digital display shows cycle parameters and time remaining.

Easy to Use

Set and save up to 10 programmable time, speed, and g-force settings using the user-friendly interface. The DASH Flex 24 arrives pre-programmed with the most common STAT settings, and can be easily reprogrammed or adjusted on the fly to meet your lab's specimen processing needs.

High Capacity STAT

The DASH Flex 24 can spin standard 75-100 mm (3-10 mL) blood tubes, 125 mm (15 mL) long tubes, and tubes with volumes up to 50 mL using a variety of different buckets - no small adapters!

Safe and Durable

Enjoy years of worry-free operation with the DASH's brushless motor and rugged, reinforced carbon fiber components. Enhance lab safety with imbalance protection and a clear, shatterproof lid, which locks until the centrifuge has stopped completely.


A Great Value

Each DASH centrifuge comes fully equipped and ready to run with a horizontal rotor, universal tube holders, and 2-year warranty included in the affordable price.


2-Stage Horizontal Rotor

Our patented swing-out horizontal rotor design incorporates a unique test tube holder that produces the ideal horizontally separated samples.

Loading Angle

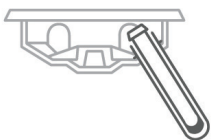

Rotation Angle

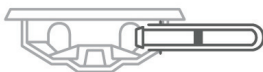

Tubes	Applications	Centrifugation	Included	Warranty
24 x Standard Tubes 12 x 125 mm/15 mL 6 x 50 mL	STAT Chemistry, Coag, PPP	Horizontal 3,450 xg Max RCF 4,500 Max RPM	Rotor, Tube Holders	2 years Rotor: lifetime warranty Extended warranty available

DASH FLEX 24



Free 30-Day Evaluation (U.S. Only)

Part Numbers	
Centrifuge	Thomas No. 20A00K823 (00-184-009-000)
Replacement Tube Holders/Buckets	02-004-0-0012K (6 qty)
50 mL Tube Holders/Buckets	02-004-0-0014K (6 qty)
Adapters, 1.5 to 2 mL	7713065 (1 qty)
1 yr Warranty Extension	03-0-0007-0003
2 yr Warranty Extension	03-0-0007-0004
3 yr Warranty Extension	03-0-0007-0005
* Custom Settings	00-100-100-101

* Add this part number to your centrifuge order to have your custom validated cycles pre-programmed at the factory

Factory Settings			
	RPM	Time	G-Force
1st Setting/Cycle	4,300	5 min	3,000 xg
2nd Setting/Cycle	4,300	6 min	3,000 xg
3rd Setting/Cycle	3,500	7 min	2,000 xg

Specifications	
Sound level	61 dBA
Dimensions (W x D x H)	14.5 x 17 x 9 in 37 x 43 x 23 cm
Weight	39 lbs (17 kg)
Power	95 to 253 VAC – 50/60 Hz



Protected by U.S. Patents #6.811.531. #D718,463S #D734,489S Other Patents Pending.



ThomasSci.com
833.544.SHIP (7447)
CustomerService@thomassci.com

Connect With Us:

